myBlend

POWER OF SYNERGIES

NUTRI + DERMA + TECH

myLEDmask*

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Please read the instructions for use carefully before use
The device is intended for cosmetic and non-medical use

myLEDmask
User's Manual

Contents

1.0	INFORMATION ABOUT THE USER MANUAL	4
2.0	SKIN AGING MECHANISMS AND INTENDED USE	4
2.1.	Skin aging mechanisms	4
2.2.	Intended use	4
3.0	LEGAL INFORMATION	4
3.1.	Legal manufacturer	4
3.2.	Liability	4
3.3.	Safety markings and symbols	5
3.4.	Precautions, safety warnings and contraindications	6
3.5.	Lifespan and warranty	7
	• Lifespan	7
	• Warranty	7
4.0	PRODUCT PRESENTATION	8
4.1.	Operating principle	8
4.2.	Description	8
	• Device components	8
	• Using the control panel	9
4.3.	Operating features	10
	• Power supply	10
	• Optical emissions delivered	10
	• Additional information	10
5.0	PRODUCT TRANSPORT	10
6.0	STORAGE	10
7.0	ASSEMBLING THE PRODUCT	10
7.1.	General information	10
7.2.	Attaching the strap	10
7.3.	Power supply	11
	• Connecting your myLEDmask to the lower part of the neckpiece	11
	• Plugging your myLEDmask into the electrical power supply	11
8.0	USING THE PRODUCT	11
9.0	FREQUENCY AND USAGE CYCLE	13
10.0	CARE AND MAINTENANCE	13
11.0	DISPOSAL	13

1.0 INFORMATION ABOUT THE USER MANUAL

This manual, which explains how to use the device, is not for sale. Reproducing it in part or in full, without permission, is prohibited. The information provided in this document may be modified without notice for technical and/or commercial reasons.

WARNING - Messages preceded by this symbol are intended to draw your attention to certain operations and situations which, if unknown to the user or incorrectly performed, may cause a risk to the user and/or to the person undergoing the facial.

2.0 SKIN AGING MECHANISMS AND INTENDED USE

2.1. Skin aging mechanisms

Aging, whether genetic, hormonal or caused by repetitive facial movements, is a natural process for humans. In the first 25 years of life, our bodies are able to fight against this phenomenon, but during the second stage of life (age 25-50), our immune system begins to weaken and generate free radicals* that harm our tissues. The functional capacities of the mitochondria** diminish, creating an imbalance between excess free radicals and insufficient antioxidant production or ingestion. This imbalance is known as oxidative stress.

We also undergo actinic*** aging caused by: the sun's UVA and UVB rays; stress factors (noise, lack of sleep, disease); the effects of climate (cold, wind, heat); air pollution (tobacco smoke, exhaust emissions, ozone, hydrocarbons, solvents); and food pollution (pesticides, detergents, pyrolyzed amino acids found in fried foods), all of which are harmful to our bodies.

All of these exogenous and endogenous factors provoke structural, functional and morphological modifications in the skin, with visible consequences including wrinkles and fine lines, slackening, puffy eyes, redness, scars and a dull complexion.

2.2. Intended use

The myLEDmask is intended to bring aesthetic benefit to the face and neck: targets wrinkles.

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Phone: +33 1 56 60 60 00

3.2. Liability

The manufacturer shall not be held liable for the safety, reliability or performance of the device if:

- The device is not used in accordance with the instructions in this user manual
- The electrical system of the place where the device is used does not meet safety standards; The manufacturer shall not be held liable for any personal injury or property damage resulting from the deterioration of the device, from a modification of the device's essential functions, from a lack of proper care or maintenance, or from the failure to comply with the instructions in this user manual.

In the event of a technical issue or other problem, contact the manufacturer:

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www.my-blend.com

Phone: +33 1 56 60 60 00

* free radicals: unstable atoms or molecules seeking to bind to other atoms, provoking chain reactions.

** mitochondria: microscopic elements, found in the cells of most living organisms, which provide energy to the cell and are essential for its survival.

*** actinic: related to or caused by light.

3.3. Safety markings and symbols

You may encounter the following safety markings and symbols in this document or on the phototherapy table.

- 

Refer to the instructions for use.
- 

This symbol is used in the user manual to warn about a potential risk for the safety of the patient or user during product use.
- 

This symbol is used to make a recommendation to the user about operating the device.
- 

Optical hazard: Infrared emissions may cause injury to the eye.
- 

Indicates that the area where this symbol is affixed is a type BF applied part.
- 

The myLEDmask by myBlend is CE-certified.
- 

Name of manufacturer
- 

Date of manufacture
- 

Item number
- 

Serial number
- IP N₁ N₁** IP protection code
- 

Direct current
- 

Play/ Pause
(The Play/ Pause symbol is used to identify the action triggered by the button: start a session or pause a session already underway.)
- 

ON/ OFF
Action of turning ON and OFF the device
ON: The device may remain ON for a maximum of 13 min 55 s
OFF: The device must remain OFF for at least one hour and a half between each session
- 

This symbol shows that pursuant to the Waste Electrical and Electronic Devices Directive (2012/19/EU) and national legislation, a product cannot be disposed of via the household waste, and must be disposed of separately within Europe.
- 

Fragile
- 

Keep away from water
- 

The charger is sensitive to electrostatic discharges. We recommend that you handle the charger with caution and ensure it is equipotential. All conductor components should be connected to earth.
- 

Temperature range: 0°C to 30°C

Solids		Water	
0	No protection	0	No protection
1	Protection from a large part of the body such as a hand and from solid objects greater than 50mm in diameter.	1	Protection from condensation and water droplets.
2	Protection against fingers or other object not greater than 80mm in length and 12mm in diameter.	2	Vertically dripping water shall have no harmful effect when the enclosure is tilted at an angle up to 15° from its normal position.
3	Protection from tools and wires greater than 2.5 millimeters.	3	Protection from water spray when at an angle up to 60° from vertical.
4	Protection from tools and wires greater than 1 mm in width.	4	Protection from water spray from any direction.
5	Protection from limited dust ingress.	5	Protection from low pressure water jets from any direction.
6	Protection from total dust ingress.	6	Protection from high pressure water jets from any direction.
Example:		7	Protection against short periods of immersion in water up to 1m depth.
		8	Protection against long, durable periods of immersion in water under pressure.

IPN₁N₂

Charger: IP 4X
myLEDmask: IP 22

Conditions for use:

- Ambient temperature: 0°C to 30°C
- Ambient humidity: 30% to 75% RH without condensation
- Pressure: 70 to 106 kPa

3.4. Precautions, safety warnings and contraindications

The mask must be unplugged and not in use during any manipulation and/or care/maintenance operations on the device.

The user must use the protective eye covers during the facial. The device is to be used while lying down.

The mask must be cleaned and disinfected after each use and after each new user.

The mask is not waterproof. Any contact with water is strictly prohibited.

Users are informed of the following risks:

In order to limit the risk of interference, portable RF (radio frequency) communication devices, including peripherals such as antenna cables and external antennae, must not be used closer than 30 cm from any part of the myLEDmask, including the cables specified by myBlend. Otherwise, the performance of the myLEDmask could be altered.

The light produced by the device does not contain any UV components and in principle presents no potential hazards. However, the large number of powerful LED lights may produce a glaring effect.

The use of accessories, transducers and cables other than those specified or provided by the device manufacturer may increase the electromagnetic emissions or decrease the resistance of the device, causing it to function inappropriately.

The near-infrared wavelength (850 nm) is not within the visible spectrum. The activity of the near-infrared LEDs therefore cannot be seen by the user. Indicator LEDs on the control panel interface show the activity of the near-infrared LEDs. Never look directly at the near-infrared LEDs when they are emitting light and, more generally speaking, do not look directly at LEDs when they are emitting light. Tests based on the NF IEC 62471 standard revealed that infrared carries a low risk (Risk Group 0) with regard to injury to the eye only.

The device must be used within a temperature range of 0° to 30° C.

Do not use the mask while it is charging. Disconnect the power supply cord before using.

The device must only be used with the material provided.

The charger must be easy to unplug.

In the event of frequent use, the LED mask must be left unused for one hour after each facial and for one hour and a half after being fully re-charged.

The charger must not be used in a humid atmosphere or near a source of water.

RISK GROUP 0

CAUTION

IR emitted from this device. Do not stare at the IR source.

Product classified to:

CE/IEC 60601-2-57-2011

Contraindications

- **Medications causing photosensitivity:** The use of retinoids or other medications causing photosensitivity is contraindicated for any phototherapy treatment. Certain fragrances and essential oils may also be contraindicated. If case of doubt, consult a dermatologist or specialist.
- **Anti-inflammatory treatments:** Reduce the benefits of LED therapy to zero.
- **Pregnancy:** The device is not to be used by pregnant women.
- **Children:** The device is not to be used by children.
- **Cancer patients:** Phototherapy can soothe many of the side effects of cancer treatments. However, the use of phototherapy on these patients should be assessed by the appropriate medical specialists.

The myLEDmask must never be covered by any material whatsoever while in use. The patient's skin must remain directly exposed to the light. As a reminder, this device is not intended for medical use, only cosmetic.

If the device enters into direct contact with the patient's intact skin, this is deemed "non-critical" (low risk of infection) and the device should be lightly disinfected. We recommend lightly disinfecting all the transparent parts of the device in contact with the skin as well as the protective eye covers before and after each facial to avoid a risk of infection otherwise.

3.5. Lifespan and warranty

• Lifespan

The estimated lifespan of the device is eight years (including average LED lifespan in the context of regular use and care), with the possible replacement of main components such as certain electronic cards.

• Warranty

The device has a two-year warranty. Our material is guaranteed for two years against any risk of dysfunction caused by a manufacturing or design defect.

The warranty is limited to the replacement of any parts recognized as defective by our technicians. We reserve the right to modify any material we have manufactured in any way we deem necessary.

The time it takes to repair the device under warranty shall not be compensated by a price reduction, by the postponement of agreed upon payment deadlines, or by an extension of the warranty.

The warranty covers the cost of labor and parts only. The warranty does not cover any dysfunction caused by the normal wear of the product, negligence, lack of care, improper use of the device or any errors in connecting or plugging in the device, including in the event of a modification to the device or intervention by technicians who have not been certified by the manufacturer.

4.0 PRODUCT PRESENTATION

4.1. Operating principle

The myLEDmask is a device that is used to perform photobiomodulation (based on light-emitting diode (LED) light) on the face and neck. Light radiates from the inside of the device onto the face and neck. This light is generated by LEDs with two different visible spectrum wavelengths: red and near-infrared. The facials are pre-programmed based on each phototype (skin tone) and recorded in the device.

4.2. Description


• Device components


The myLEDmask is composed of a LED mask and a battery charger. A **control panel located on the outside of the neckpiece** is used to operate the device.

The battery is also housed in the neckpiece. When fully charged, it enables the device to be used several times.

The myLEDmask features **plastic eye protections on the inside** which protect your eyes during use. These protections also enable you to correctly position the mask on your face. For each eye there is also a cover piece which can be removed and cleaned.


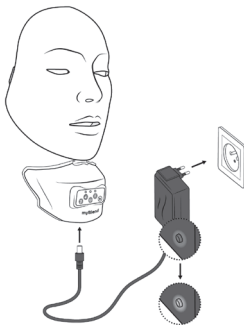
The retaining strap allows you to attach the mask around your head.

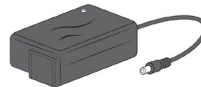
Charger

The charger cord connects to the lower part of the neckpiece. (see 7.3 Power supply).



When plugged into an electrical network, **the mask charges in 3.5 hours** (cord length: 1800 mm). Use one of the **four specific adapters** provided, based on your country.

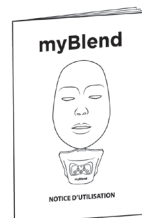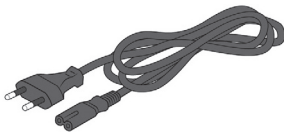The User Manual contains all the necessary legal information, product information and instructions for use.

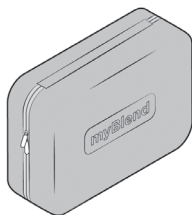
The storage case for the myLEDmask contains all the accessories and the user manual, for easy transportation.

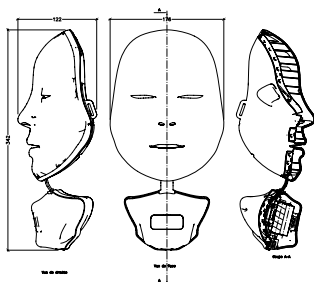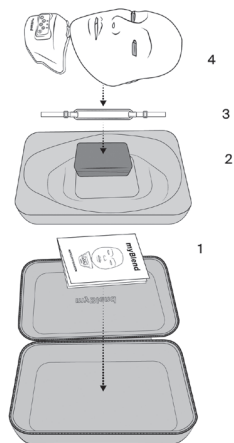myLEDmask characteristics:


Size: 17,6 x 34,2 x 12,2 cm
Weight: 650 g

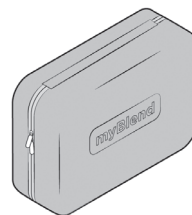Storage case characteristics:

Size: 36,5 x 22,5 x 15cm
Weight: 623 g

• Using the control panel

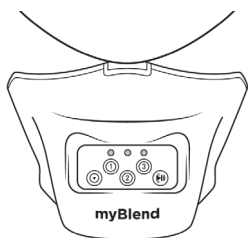Turning on the device

Press the button that corresponds to your phototype to launch your "regenerating experience" program.

- Phototype 1: light, pale white, white, fair
- Phototype 2: medium, white to olive, moderate brown
- Phototype 3: brown, dark brown, very dark, black

Start/ Stop and Resume program

Your myLEDmask will automatically stop at the end of the session (see 8. Using the product)

4.3. Operating features

- **Power supply – Charger: MASCOT charger model number 2116. The use of any other charger is strictly prohibited.**
 - Voltage: 100 to 240 volts, AC
 - Frequency: 50-60Hz
 - Intensity: 0.35 A

The charger is sensitive to electrostatic discharges. We recommend that you handle the charger with caution and ensure it is equipotential. All conductor components should be connected to earth.

- **OPTICAL EMISSIONS DELIVERED**
 - Red LEDs= 630 nm
 - Number of red LEDs on face: 122
 - Number of red LEDs on neck: 22
 - Infrared LEDs= 850 nm
 - Number of infrared LEDs on face: 122
 - Number of infrared LEDs on neck: 22

- **Precise operating characteristics**

Explanation of the 3 defined programs

Phototype	Duration	Battery autonomy	Number of complete programs
Program 1 Phototype 1 (light, pale white, white, fair)	5 min 35	35 min	6
Program 2 Phototype 2 (medium, white to olive, moderate brown)	11 min 10	35 min	3
Program 3 Phototype 3 (brown, dark brown, very dark, black)	13 min 55	35 min	2

5.0 PRODUCT TRANSPORT

Atmospheric conditions: $0^{\circ}\text{C} \leq 0 \leq +45^{\circ}\text{C}$. The myLEDmask comes in a storage case specially designed for safe transportation.

6.0 STORAGE

The device must be stored in a dry place with no direct sunlight. During storage, humidity must be between 10% and 80% with no condensation, and the temperature between 0°C and 35°C , with atmospheric pressure ranging from 70 kPa to 106 kPa. Make sure your myLEDmask is stored in its case when not in use. This will protect it from exposure to light and from potentially damaging impacts or falls.

7.0 ASSEMBLING THE PRODUCT

7.1. General information

The myLEDmask device is delivered in one piece, already assembled.

Before using it for the first time, it is recommended that you fully charge the mask for at least **3.5 hours using the charger**. When the battery is fully charged, the indicator light on the charger turns green, meaning you must unplug the charger. **Let the mask sit for one hour and a half after each complete charging**, and ensure it is unplugged before using.


The device must be unplugged when in use.


To connect the device to the power supply and turn it on, an IEC electrical outlet with two terminals is all that is needed (see 4.2 Device components).

Device settings are adjusted ahead of delivery. There is therefore no need for any calibrations by the user. However, it is important to follow the instructions in **Chapter 8 – Using the Product**.


7.2. Attaching the strap

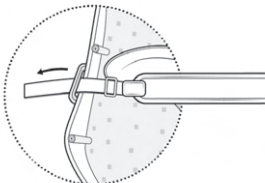Introduce the strap into the slotted tab on either side of the mask.

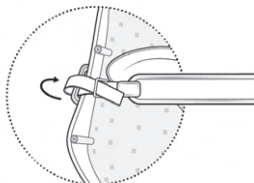
Slide the strap into the plastic loop.

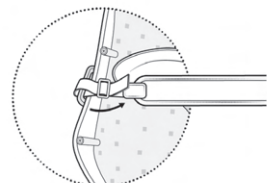Adjust the strap so it fits snugly around your head.

7.3. Power supply

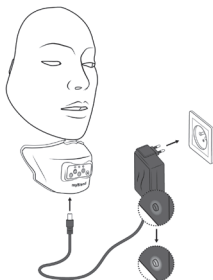• Connecting your myLEDmask to the lower part of the neckpiece

To operate your myLEDmask, you will need to connect to a 100-240V 50Hz electrical power supply.


The charger must not be used in a humid atmosphere or near a source of water.

Connect the charger cord's jack plug to the lower part of the neckpiece.

• Plugging your myLEDmask into the electrical power supply

Plug the charger into an electrical outlet (100-240V 50Hz) using the specific adaptor suited to your country (four different adaptors are provided).

Remove the charger cord from the neckpiece and from the electrical outlet, and let the mask sit for at least one hour and a half after each full charging.

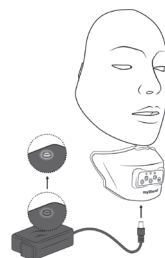Check the state of charge of the LED mask.

When you plug in your LED mask, an indicator light shows the battery's state of charge.

Indicator light	State
Orange	Charging
Green	Fully charged

Note:

Unplug the mask when the green light appears, and let it sit for at least one hour and a half before using.

Your fully charged myLEDmask can then be used for up to six sessions, with one hour pause between each session (see 4.3 Operating features). To accommodate regular use thereafter, simply recharge the device as needed using the charger cord. Once your myLEDmask is charged, you must remove the charger cord from the neck piece before using the device.

The device must be unplugged when in use.

8.0 USING THE PRODUCT

No particular skills are needed to use the myLEDmask. It was designed for ease of use thanks to predefined programs, which minimize the number of steps required.

STEP 1: Prepare your skin for the facial

Clean and dry your face before each session in order to eliminate makeup and impurities.

Refer to step 8 to learn how to clean your myLEDmask

STEP 2: Turn on your myLEDmask


Turn on your mask by pressing the "On" button  found on the control panel in the mask neckpiece.

The light above Button 1 on the control panel will come on to confirm that the device has been turned on. Phototype 1 is automatically preselected.


STEP 3: Select your phototype

Choose your facial program. To select your phototype*, use the following guidelines:


Phototype I:
light, pale white, white, fair.
Left-hand button.
This prototype is selected by default.

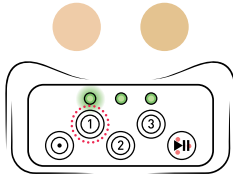Phototype II:
medium, white to olive, moderate brown.
Middle button.

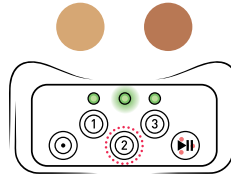Phototype III:
brown, dark brown, very dark, black.
Right-hand button.

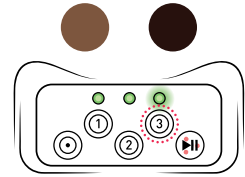*The facial is based on your phototype. Choosing a facial that does not match your phototype will not reap the full benefits of this device.

STEP 4: Position the mask on your face

Place the mask on your face and adjust the strap to ensure the mask fits snugly.

STEP 5: Lie down

To fully enjoy the benefits of your facial, we recommend that you lie on your back, keep your eyes closed and adjust the eye protections to avoid any glare.

STEP 6: Begin your facial

Once the phototype has been selected (or preselected if using Phototype I), press Start (raised button) for two seconds:

 Start/ Stop and Resume to begin your facial

During the session, the light of the phototype selected will blink to confirm the choice of phototype.


STEP 7: End your myLEDmask facial/ Turn off the device

Your LED mask will automatically turn off at the end of your facial.

You may interrupt the session at any time by pressing the Start/ Stop & Resume program button  for two seconds.

If you pause for more than 30 seconds by pressing the Start/ Stop & Resume program button , the indicator light will blink slowly for 30 seconds and the mask will automatically be reprogrammed to Phototype I.

If you pause for more than 60 seconds, the device will automatically turn off for safety reasons. It is preferable not to stop a facial in progress.

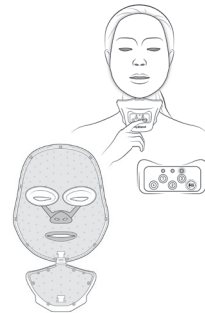STEP 8: Cleaning & disinfecting myLEDmask (using 70% IPA)

At the end of each session, clean your myLEDmask and place it in its storage case with the charger.

Clean the mask's inside face and neck surfaces that are in contact with the skin using a fresh individual hand wipe or an alcohol-based solution (70% concentration).

DO NOT USE ANY other cleaning products, which might damage the plastic.

9.0 FREQUENCY AND USAGE CYCLE

A complete cycle consists in several sessions spaced at least 48 hours apart (3 to 4 sessions a week), and lasts for six to eight consecutive weeks. This cycle may be repeated with a one-month pause between each cycle.

10.0 CARE & MAINTENANCE

Using the mask for any purpose other than that described in this manual is prohibited and will result in the automatic cancellation of the warranty. Should any problems arise, contact our After-Sales Service or your retailer. Any modification to the device will automatically render the warranty null and void.

11.0 DISPOSAL

Disposal concerns either the entire product at the end of its lifecycle, or parts that have become defective or damaged over time. At the end of its lifecycle, the device must be disposed of in compliance with the regulations regarding waste management in force in the country in question. Feel free to contact myBlend if you have any questions about the applicable regulations. The device contains plastic, mechanical and electronic components as well as a battery.

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MYBLEND SAS