myBlend

myLIFTINGroller

Notice d'utilisation
Instruction manual
Benutzerhandbuch
Gebruikershandleiding
Manuale utente

FRANÇAIS

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ENGLISH

COMPANY INFORMATION

myBlend is committed to revolutionizing skincare through innovation, sensoriality and beauty-tech. Based on the principles of epigenetics, our mission is to enable every woman to reveal her skin's true potential.

GENERAL DEVICE INFORMATION

DESCRIPTION OF THE DEVICE

Thanks to its advanced technology, myLIFTINGroller fights skin slackening by stimulating facial muscles, boosting skin firmness and elasticity. Through in-depth action, it helps redefine facial contours while improving skin quality and natural radiance.


It combines 3 mechanisms: mechanical massage, electrostimulation and vibration for advanced results. Features are enhanced, cheekbones lifted.

Integrated into a daily routine, this powerful, intuitive tool transforms every gesture into a comprehensive treatment, tailored to your skin's needs.


DEVICE SPECIFICATIONS

Operating temperature	0°C to 40°C 32°F to 104°F
Storing and transport temperature	-5°C to 55°C 23°F to 131°F
Operating humidity	30% to 95% without condensation
Storage and Transport Humidity	20% to 95% without condensation
Charging supply	5V 2A
Charging port	USB-C
Charging cable	USB-C to USB-A
Charging cable length	1 meter / 3.28 feet
Battery	3.7V, 400 mAh
Battery type	Lithium-ion
Charging time	120 minutes / 2 hours
Treatment time	7 minutes
Dimensions of the device	157.9 x 80.1 x 37.3 mm
Weight of the device	140g / 4.94oz


DEVICE ACCESSORIES

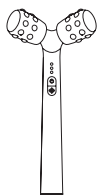
Device

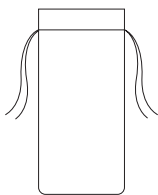Travel
Pouch

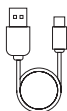Charging
cable

User
Manual

Quick Start
Guide

ADVERSE REACTION

If you have any adverse reaction, such as discomfort, skin discoloration, skin irritation or other concerns, from using this Device, stop using it immediately and consult your Healthcare Provider.

Skin irritation, bruising, contact burns, dizziness, nausea, and headaches are rare but possible side effects. Immediately stop using the device and contact your Healthcare Provider.

Skin irritation and burns beneath the electrodes have been reported with the use of other powered muscle stimulators. If you experience skin irritation or burning, discontinue use.

If you experience a persistent headache and/or other persistent painful sensation during use of this device, stop using the device and consult with your Healthcare Provider.

CONTRAINDICATIONS

- DO NOT use the device if you have any active cancer or cancerous lesion.
- DO NOT use the device if you are pregnant, planning on becoming pregnant, or breastfeeding.
- DO NOT use the device on people under the age of 18.
- DO NOT use the device if you have severe acne.
- DO NOT use the device over areas with phlebitis, thrombophlebitis, and varicose veins.
- DO NOT use the device if you have a known allergic reaction to zinc alloy. The device contains plastics and metals that come into contact with the skin during use. When in doubt of allergies first contact your Healthcare Provider before using the device.
- DO NOT use the device if you have an electronic implanted device such as a defibrillator neurostimulator, pacemaker or ECG monitor. If you use a body-worn device such as an insulin pump, consult your Healthcare Provider.
- DO NOT use the device on the thyroid gland in the neck area, as it may stimulate its growth.
- DO NOT use the device on the neck if you have a known history of carotid artery disease, stroke or transient ischemic attack (TIA), carotid stenosis, unstable blood

pressure, easy fainting, thyroid disease or prior to surgical procedures in the neck area.

Note: We highly suggest taking caution when taking any (skin related) damaging medications. If in doubt, consult your Healthcare Provider.

Note: If you are under the care of a Healthcare Provider, consult your Healthcare Provider before using this device.

WARNINGS

- STOP using the device and consult your Healthcare Provider if you experience an adverse reaction from using the device.
- STOP using the device if you experience any discomfort or difficulty breathing or swallowing occurring from application of the device on the neck.
- STOP using the device if you experience symptoms of hypotension (dizziness or light-headedness) while using the device.
- STOP using the device if it is overheating. Turn off the device immediately.
- STOP your treatment if it feels hot or your skin feels uncomfortably warm.
- DO NOT let children play with the device. Cleaning and charging the device shall not be done by children.
- DO NOT use the device to treat rosacea, moles, warts, open sores, cancerous lesions, or any skin condition.
- DO NOT use the device on spots without pigmentation or areas on the skin that have vitiligo.
- DO NOT place the device if you have any open wounds.
- DO NOT use the device outdoors to prevent damage to the product.
- DO NOT sit, lay down or step on the device. This will break the device.
- DO NOT drop or insert any object into any openings of the device.
- DO NOT use the device if it is not working properly, or if it has been dropped or damaged. For warranty and repairs contact customer support.
- DO NOT use or store the device near hot surfaces.
- DO NOT charge the device when it is close to or might get in contact with water.
- DO NOT submerge the device with water.
- DO NOT store the device in extreme heat or cold.
- DO NOT place the device in direct sunlight.
- DO NOT attempt to modify the device in any way. Modifying or opening the device will void your warranty. the device contains no serviceable parts. This device contains batteries that are non-replaceable
- DO NOT perform maintenance or take apart the product. In case of malfunctioning or damage, contact customer support.
- DO NOT use the device if damaged or when the charging cable is damaged.
- DO NOT use the device with non-manufacturer-approved accessories.
- DO NOT use solvents or abrasive cleaners on the device.
- DO NOT use the device in any other way then specified in this manual.
- DO NOT use the device if you have had a recent surgery or other surgical procedures.
- DO NOT use the device over a bacterial infection as this may result in the worsening of symptoms.
- DO NOT use the device on painful areas. Consult your Healthcare Provider when you experience painful area's during use.
- DO NOT use the device on, near, or by children or those with reduced physical and mental abilities. Close supervision is necessary when this device is used, cleaned or stored near children or those with reduced physical or mental abilities.
- DO NOT use the device with any other USB-C cable than provided by myBlend.
- DO NOT use the device with any other topical other than specified by myBlend.

- DO NOT use the device over areas of skin that lack normal sensation. Only use the device on normal, intact, clean and healthy skin.
- DO NOT use the device while driving, operating machines or during an activity in which voluntary muscle contractions may put users or the public to any risk of injury.
- DO NOT use the device on swollen, infected or inflamed area's or skin eruptions.
- DO NOT use the device in the eye area (circular muscle within the orbital rim), on the temples, the breast or groin area.
- DO NOT apply excessive pressure when using the device on the anterior neck.
- DO NOT use the device directly on the midline of the neck, specifically avoiding the thyroid.
- DO NOT apply stimulation across your chest as the introduction of electrical current into the chest may cause rhythmic disturbances to your heart, which could be lethal.
- DO NOT use the device on piercings or other metal implants
- DO NOT use the device at a place where other devices could interfere by emitting electromagnetic energy. This device can be negatively affected by interference of other devices.
- DO NOT operate the device in the presence of electronic monitoring equipment (e.g. heart monitors, heart tracing or close to (1m) shortwave devices like micro-waves and electromagnetic therapy equipment).

Note: The device should only be supplied with 5V DC low voltage. If any part of the main adapter or cable is damaged, discontinue use. Contact myBlend to arrange a cable replacement.

PRECAUTIONS

- Do not scratch or damage the metal node surfaces.
- It will void your warranty if you do not operate the device in accordance with the instruction in the manual.
- Read all instructions carefully before using this device.
- Store the device according to the protocol worked out in this manual.
- This device is intended to be used by one person; we recommend to not share it with anyone else for reasons of hygiene.
- The device is designed for cosmetic use only, for individuals with good health. Please contact your Healthcare Provider if you are unsure about your health status.
- The battery of the device is not removable by the user. If the device is to be disposed of, do not put into a fire. Dispose of the device per local regulations.
- Inspect the device and do not use the device if you see any signs of damage such as internal battery leaks. Contact our customer support.
- The device should only be supplied with 5V DC low voltage. If any part of the main adapter or cord is damaged, discontinue use. Contact customer support to arrange a replacement.
- Consult your Healthcare Provider before using the device if you have a fever, acute inflammation, a severe circulatory disorder, anaesthetic, or desensitized areas.
- We highly suggest taking caution when taking any medications. If in doubt, consult your Healthcare Provider.
- The long-term effects of transcutaneous electrical stimulators are unknown. If you have suspected or diagnosed heart disease, you should follow precautions recommended by your Healthcare Provider.
- Some users may experience skin irritation or hypersensitivity due to the electrical stimulation or serum.
- Skin irritation and burns beneath the electrodes have been reported with the use

of other powered muscle stimulators. If you experience skin irritation or burning, discontinue use immediately.

- If using the device causes persistent pain or discomfort, discontinue use immediately.
- A slight tingling sensation may occur while using the device. Lowering the intensity may reduce or eliminate this sensation.
- There may be a perception of flashing lights during the use of the device, resulting from the stimulation of your optic nerve. If you notice this condition constantly without the use of your device, your Healthcare Provider.
- Some users may experience a metallic taste in their mouth during the use of your device, resulting from stimulation and varying sensitivity around the mouth. If you notice this condition constantly without the use of your device, consult your Healthcare Provider.

TRAVEL AND INTERNATIONAL USE

The device's lightweight design makes it a perfect travel companion. Here are some travel tips for you:

- If you are travelling to a different region of where you bought the Device, please assure that the device is only supplied with a 5V DC Low Voltage.
- Confirm with the steward that the Device can be carried and/or used on the airplane prior to flight.

DISCLAIMER

This Device is not intended to diagnose, cure, or prevent any disease. myBlend makes no claims, representations, or warranties regarding the ability of this Device to cure any physical, skin, or mental condition. A qualified Healthcare Provider should be consulted regarding any condition requiring medical attention.

POTENTIAL RISKS

Immediately after treatment skin may appear red. Skin should return to normal within 1-2 hours.


Temporary headaches immediately following the treatment. Headaches can be treated with an appropriate pain killer and rehydration.

HOW TO USE THE DEVICE


GENERAL USE

1. Prior to using the device, make sure it is fully charged.
2. Clean the skin that you want to treat.
3. Apply myBlend Firming Superserum or a water-based serum to drive the current through your skin.
4. Make sure the device is unplugged and place the device on the skin.
5. Press and hold the ON/OFF to turn on the device, first LED will turn on.
 - 5.1. First Microcurrent LED will turn on. In case of low battery, the light next to the charging port will start blinking.


- 5.2. To set the desired microcurrent level, press the ON/OFF button briefly. Each press cycles through the available levels.
- 5.3. Gently glide the device over the treatment area in slow, upward or outward motions, ensuring consistent contact with the skin.
- 5.4. To control the motor vibration feature, press the vibration button briefly. Each press activates or deactivates will deactivate the vibration, allowing you to adjust the session to your preferences.
- 5.5. Keep the device constantly moving, without holding it in the same spot for several seconds.
- 5.6. The device switches off after a 7-minute treatment automatically. You may stop the treatment at any time by pressing the ON/OFF button.

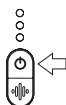
5.0

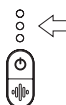5.1

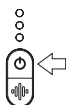5.2

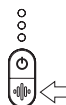5.4

6. Take the device off the treatment area.
7. Clean the Device using a water-based wipe or a damp clean cloth. Ensure the Device is NOT connected to the power supply when cleaning.
8. Repeat these steps based upon the treatment protocol described in the next part.

TREATMENT PROTOCOL

Use the device 3-4 times per week, based upon your lifestyle and skin condition. It is not recommended to use the device more than 7 days per week. Do not use the device more than once per day.

CHARGING

Charge myLIFTINGroller fully for 2 hours before your first use. Only use a 5V adapter or computer USB port (avoid fast chargers as they can damage the battery).

To charge the device insert the USB-C cable into the device. Insert the USB-A side of the cable into a USB power supply. The LED indicator on the device will start to flash slow or blink slow to indicate charging. When the device is fully charged the LED indicator will turn off.

Before charging, make sure that the plug and socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.

Note: Always unplug the device before use, the device will not work while being charged.

CARE AND STORAGE

- The lifetime of the Lithium battery can be maximized by allowing it to fully charge and then fully discharge during regular use.
- Store the device at room temperature in dry conditions.
- Store the device in the supplied case or the packaging. Make sure the device will

- not bend, stretch or crush while being stored or when travelled with.
- Avoid storing the device in dusty conditions, in direct sunlight or in damp and humid conditions.
 - This device contains no serviceable parts.
 - Before each use, with the device switched off, clean the Device following the cleaning instructions below.

IMPORTANT

- The product has a limited lifecycle of two years under nominal use. Following the treatment protocol for two years is considered as nominal use.
- The device must be handled with care.

CLEANING

Clean the Device with a water-based wipe or damp clean cloth. Stubborn stains can be removed by wiping gently with 70% alcohol.

- DO NOT use any other cleansing Devices as these may damage the Device.
- DO NOT clean the device while it is turned on.
- DO NOT clean the device while charging.

TROUBLESHOOTING

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please return to myBlend.

- DO NOT attempt to investigate any further.
- DO NOT attempt to modify the device in any way. The device has NO parts you can repair yourself.

Problem	Cause	Solution
The device doesn't switch on.	The device is not fully charged, because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The device is not fully charged, because the charging cable was not properly inserted into the device.	Insert the charging cable properly into the device and charge the device.
	The power outlet didn't work during charging.	Try a different power outlet.
	The on/off button has not been pushed one full second.	Push the on/off button for 1 second straight.
The device switches off unexpectedly.	The battery level is too low.	Charge the device.

DISPOSAL

DO NOT dispose of the device with household waste. The Device should be recycled. The Device contains a Lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.

WARRANTY

The device is warranted to be free from defects in material and workmanship for 24 months from its original purchase date under normal use. This warranty extends only to the original retail purchaser with original proof of purchase and only when purchased from an authorized retailer or reseller. This warranty only extends if the device is used in conjunction with authorized Device accessories. If the Device should become defective within the warranty period, contact the retailer that you purchased the device. This warranty does not cover Devices damaged by:

- Accident, misuse, abuse, or alteration.
- Use with unauthorized accessories or use other than as instructed.
- Connecting to incorrect current and voltage.

MyBlend shall not be responsible for any incidental, special or consequential damages resulting from the use of this Device. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to two years from date of original purchase.

This warranty gives you specific legal rights, and you may also have other rights that vary from country to country. Some countries do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you.

This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose are hereby disclaimed.

If the Device is defective within the warranty period please contact the Customer Service Team.

To register your Device and be kept up to date with new and exciting developments from the Customer Service team, please take a moment to sign up on our website www.my-blend.com. To complete your registration online you will need to enter the Device batch number.

SAFETY

The Device has been tested to or complies with the following standards:

IEC/EC 62133:2017 2nd ed. Secondary cells and batteries containing alkaline or other non-acid electrolytes - Safety requirements for portable sealed secondary cells, and for batteries made from them, for use in portable applications - Part 2: Lithium systems

EN 55014-1:2006+A2:2011 Electromagnetic compatibility. Requirements for household appliances, electric tools and similar apparatus. Emission

EN 55014-2:2015 Electromagnetic compatibility. Requirements for household appliances, electric tools and similar apparatus. Immunity. Device family standard. Restriction of Hazardous Substances (RoHS) Directive (2011/65/EU)

Registration, Evaluation, Authorisation and Restriction of chemicals (REACH) directive (EC 1907/2006)

IEC 60335-1:2020 – Household and similar electrical appliances - Safety - Part 1: General requirements

EN IEC 60335-1:2023 + A11:2023 – Household and similar electrical appliances – Safety – Part 1: General requirements (European Standard incorporating Amendment A11)

EN 62233:2008 – Measurement methods for electromagnetic fields of household appliances and similar apparatus with regard to human exposure

FAQ

For more information, visit the FAQs section in our website.

SYMBOLS

Indicates the Device is CE Marked and authorized to sell in Europe. (Lower-Voltage-Directive)

The Waste of Electrical and Electronic Equipment (WEEE), Indicates the Device may not be discarded as household waste (2012/19/EU).

Restriction of hazardous substances.

This warning symbol on the device label is a warning to read the enclosed documents before use and take all safety precautions required.

This symbol on the device label means: Read the user manual before you start using the device.

This symbol means made in China.

This symbol means the upper and lower operating temperatures that the device can be safely used at.

This symbol means the upper and lower operating humidity that the device can be safely used at.

This symbol means that the cardboard packaging that the device came in can be recycled.

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Light Tree Ventures Europe B.V.