myBlend

myLEDmask 2

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Instruction manual

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Irish Gaelic

English	4
Français	15
Deutsch	26
Nederlands	37
Italiano	48
Español	59
Portuguese	70
Česky	81
Polski	92
Gaeilge	103

ENGLISH

THE ULTIMATE BEAUTY TECH EXPERIENCE

Enhance your skin with the new generation myLEDmask 2, a perfect combination of **288 red and infrared LEDs** for complete and targeted action on the skin.

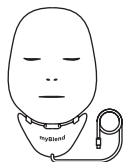
The result of **clinical research** and **consumer studies**, myLEDmask 2 transcends the beauty routine and offers an incomparable wellness experience.

Exclusive design, cutting-edge technology: **excellence at the service of your beauty.**

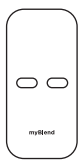
SPECIFICATIONS OF myLEDmask 2

LED wavelengths	630nm±5nm (Red), 850±5nm (Infrared)
Max. output optical power density	Red: 16mW/CM ² , Infrared: 7.5mW/CM ²
LEDs	144 Dual-core LEDs (288 diodes)
Operating temperature	0°C to 30°C 32°F to 86°F
Storing and transport temperature	-5°C to 55°C 23°F to 131°F
Operating humidity	30% to 95% without condensation
Storage and Transport Humidity	20% to 95% without condensation
Charging supply	5V 2A
Charging port	USB-C
Charging cable	USB-C to USB-C
Charging cable length	1 meter / 3.28 feet
Battery	1500mAH (2PCS)
Battery type	Lithium-ion
Charging time	120 minutes / 2 hours
Recommended use	3 to 4 times per week for 8 weeks, keep at least 24 hours between sessions
Treatment time	Mode 1: 5:30 min Mode 2: 11 min Mode 3: 14 min
Dimensions of the device	myLEDmask 2: 297 x 177 x 120,5 mm Controller: 133,5 x 64 x 21 mm
Weight of the device	myLEDmask 2: 610 grams Controller: 150 grams

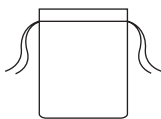
DEVICE ACCESSORIES

myLEDmask 2

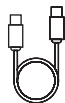
Controller

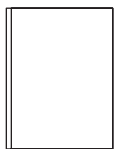
Pouch


Travel case

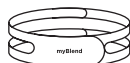
Charging cable

Instruction manual

Quick Start Guide

Strap

Eye protection

ADVERSE REACTION

If you have any adverse reaction, such as discomfort, skin discoloration, skin irritation or other concerns, from using this Device, stop using it immediately and consult your Healthcare Provider.

Skin irritation, bruising, contact burns, dizziness, nausea, and headaches are rare but possible side effects. Immediately stop using the device and contact your Healthcare Provider.

CONTRAINDICATIONS

- DO NOT use the device if you have any active cancer or cancerous lesion.
- DO NOT use the device if you are pregnant, planning on becoming pregnant, or breastfeeding.
- DO NOT use the device on people under the age of 18.
- DO NOT use the device if you have severe acne.
- DO NOT use the device over areas with phlebitis, thrombophlebitis, and varicose veins.
- DO NOT use the device if you have a known allergic reaction to plastics. The device contains plastics that come into contact with the skin during use. When in doubt of allergies first contact your Healthcare Provider before using the device.
- DO NOT use the device if you suffer from any photosensitive disorder (sensitization to light).
- DO NOT use the device if you suffer from any skin condition.
- DO NOT use the device if you are taking any medication or drugs that can cause photosensitivity. If you are not sure about your medication contact your Healthcare Provider.
- DO NOT use the device if you suffer from light-induced headaches
- DO NOT use the device if you have a medical history of seizures or epilepsy that was triggered by light.

Note: We highly suggest taking caution when taking any (skin related) damaging medications. If in doubt, consult your Healthcare Provider.

Note: If you are under the care of a Healthcare Provider, consult your Healthcare Provider before using this device.

WARNINGS

- STOP using the device and consult your Healthcare Provider if you experience an adverse reaction from using the device.
- STOP using the device if you experience any discomfort or difficulty breathing or swallowing occurring from application of the device on the neck.
- STOP using the device if you experience symptoms of hypotension (dizziness or light-headedness) while using the device.
- STOP using the device if it is overheating. Turn off the device immediately.
- STOP your treatment if it feels hot or your skin feels uncomfortably warm.
- DO NOT let children play with the device. Cleaning and charging the device shall not be done

by children.

- DO NOT use the device to treat rosacea, moles, warts, open sores, cancerous lesions, or any skin condition.
- DO NOT use the device on spots without pigmentation or areas on the skin that have vitiligo.
- DO NOT shine the light from the LEDs directly into your eyes or someone else's eyes. This may cause temporary loss of sight.
- DO NOT place the device if you have any open wounds.
- DO NOT use the device outdoors to prevent damage to the product.
- DO NOT sit, lay down or step on the device. This will break the device.
- DO NOT drop or insert any object into any openings of the device.
- DO NOT use the device if it is not working properly, or if it has been dropped or damaged. For warranty and repairs contact customer support.
- DO NOT use or store the device near hot surfaces.
- DO NOT charge the device when it is close to or might get in contact with water.
- DO NOT submerge the device with water.
- DO NOT store the device in extreme heat or cold.
- DO NOT place the device in direct sunlight.
- DO NOT attempt to modify the device in any way. Modifying or opening the device will void your warranty. The device contains no serviceable parts. This device contains batteries that are non-replaceable.
- DO NOT perform maintenance or take apart the product. In case of malfunctioning or damage, contact customer support.
- DO NOT use the device if damaged or when the charging cable is damaged.
- DO NOT use the device with non-manufacturer-approved accessories.
- DO NOT use solvents or abrasive cleaners on the device.
- DO NOT use the device in any other way than specified in this manual.
- DO NOT use the device if you have had a recent surgery or other surgical procedures.
- DO NOT use the device over a bacterial infection as this may result in the worsening of symptoms.
- DO NOT use the device on painful areas. Consult your Healthcare Provider when you experience painful areas during use.
- DO NOT use the device on, near, or by children or those with reduced physical and mental abilities. Close supervision is necessary when this device is used, cleaned or stored near children or those with reduced physical or mental abilities.
- DO NOT use the device with any other USB-C cable than provided by myBlend.
- DO NOT use the device with any other topical other than specified by myBlend.
- DO NOT use the device over areas of skin that lack normal sensation. Only use the device on normal, intact, clean and healthy skin.
- DO NOT use the device on swollen, infected or inflamed areas or skin eruptions.
- DO NOT apply excessive pressure when using the device on the anterior neck.

Note: The device should only be supplied with 5V DC low voltage. If any part of the main adapter or cable is damaged, discontinue use. Contact myBlend customer support (serviceclient@my-blend.com) to arrange a cable replacement.

PRECAUTIONS

- Do not scratch or damage the light-emitting surface.
- It will void your warranty if you do not operate the device in accordance with the instruction in the manual.
- Read all instructions carefully before using this device.
- Store the device according to the protocol worked out in this manual.
- This device is intended to be used by one person; we recommend to not share it with anyone else for reasons of hygiene.
- The device is designed for cosmetic use only, for individuals with good health. Please contact your Healthcare Provider if you are unsure about your health status.
- The battery of the device is not removable by the user. If the device is to be disposed of, do not put into a fire. Dispose of the device per local regulations.
- Inspect the device and do not use the device if you see any signs of damage such as internal battery leaks. Contact our customer support.
- The device should only be supplied with 5V DC low voltage. If any part of the main adapter or cord is damaged, discontinue use. Contact customer support to arrange a replacement.
- Consult your Healthcare Provider before using the device if you have a fever, acute inflammation, a severe circulatory disorder, light sensitivity, anaesthetic, or desensitized areas.
- We highly suggest taking caution when taking any medications. If in doubt, consult your

Healthcare Provider.

TRAVEL AND INTERNATIONAL USE

The device's lightweight design makes it a perfect travel companion. Here are some travel tips for you:

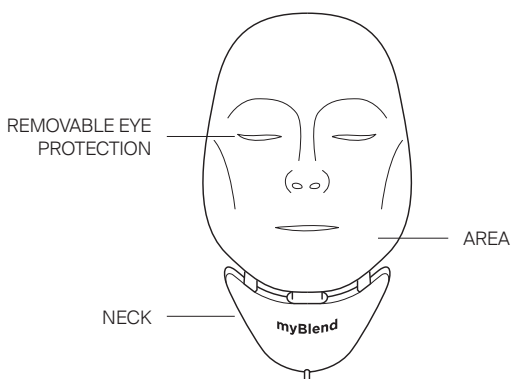
- If you are travelling to a different region of where you bought the Device, please assure that the device is only supplied with a 5V DC Low Voltage.
- Confirm with the steward that the Device can be carried and/or used on the airplane prior to flight.

DISCLAIMER

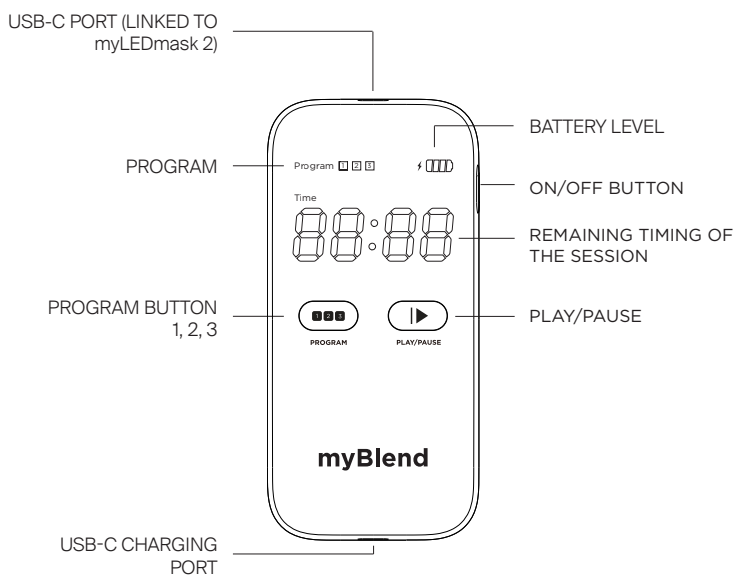
This Device is not intended to diagnose, cure, or prevent any disease. myBlend makes no claims, representations, or warranties regarding the ability of this Device to cure any physical, skin, or mental condition. A qualified Healthcare Provider should be consulted regarding a condition requiring medical attention.

PRODUCT OVERVIEW

MASK

CONTROLLER

HOW TO USE myLEDmask 2

LIGHT SENSITIVITY TEST

We highly recommend users perform a light sensitivity test before their first use of the device.

1. Place the device on the inner forearm
2. Turn on the device and start a treatment on 'Mode 1'
3. Wait till the treatment is done. The treatment will take 5 minutes and 30 seconds to complete.
4. If your skin feels hot or uncomfortable during the treatment directly stop! You might be sensitive to light.
5. The device will automatically turn off after the treatment is done. Remove the device and wait for at least 24 hours.

Check your skin regularly over those 24 hours, if your skin appears red and/or itchy the device is unsuitable for you. In the case of light sensitivity please contact customer support.

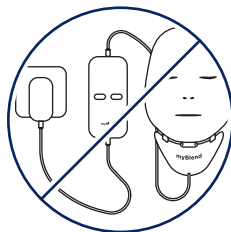
CHARGING

Charge your myBlend Mask fully for 2 hours before your first use. Only use the accompanied power adapter when charging the device.

CHARGER  To charge the device insert the USB-C cable into the controller.

The battery indicator on the controller will start cycling through the charging levels to indicate charging.

When the device is fully charged, the battery icon will display full charge.

Before charging, make sure that the plug and socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.

Note: Always unplug the device before use, the device will not work while being charged.

GENERAL USE

1. Prior to using the device, make sure it is fully charged.
2. Clean the skin that you want to treat.
3. Make sure the controller is unplugged from the main voltage and attach the cable of the mask in the controller. Place the device on the face. Secure the mask with the use of the accompanied strap.
4. Press and hold the ON/OFF button on the controller to turn on the device.
 - a. To select the desired program, click the 'program' button on the controller to cycle through the modes.
 - i. Mode 1: 5:30 minutes
 - ii. Mode 2: 11 minutes
 - iii. Mode 3: 14 minutes

- b. Press the 'play' button to start your session. The time indications will be displayed on the controller.
 - c. The device switches off after the treatment time is completed. You may pause the treatment by pressing the 'play' button on the controller.
5. Take the device off your face.
 6. Clean the Device using a water-based wipe or a damp, clean cloth. Ensure the Device is NOT connected to the power supply when cleaning.
 7. Repeat these steps based on the treatment protocol described in the next part.

CARE AND STORAGE

- The lifetime of the Lithium battery can be maximized by allowing it to fully charge and then fully discharge during regular use.
- Store the device at room temperature in dry conditions.
- Store the device in the supplied case or the packaging. Make sure the device will not bend, stretch or crush while being stored or when travelled with.
- Avoid storing the device in dusty conditions, in direct sunlight or in damp and humid conditions.
- This device contains no serviceable parts.
- Before each use, with the device switched off, clean the Device following the cleaning instructions below.

CLEANING

- Clean the Device with a water-based wipe or damp clean cloth. Stubborn stains can be removed by wiping gently with 70% alcohol.
- DO NOT use any other cleansing Devices as these may damage the Device.
- DO NOT clean the device while it is turned on.
- DO NOT clean the device while charging.

TROUBLESHOOTING

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please return to myBlend.

- DO NOT attempt to investigate any further.
- DO NOT attempt to modify the device in any way. The device has NO parts you can repair yourself.

Problem	Cause	Solution
The device doesn't switch on.	The device is not fully charged, because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The device is not fully charged, because the charging cable was not properly inserted into the device.	Insert the charging cable properly into the device and charge the device.
	The power outlet didn't work during charging.	Try a different power outlet.
	The on/off button has not been pushed one full second.	Push the on/off button for 1 second straight.
	The cable of the mask is not properly connected in the controller.	Insert the cable of the mask fully in the controller.
The device switches off unexpectedly.	The battery level is too low.	Charge the device.

DISPOSAL

DO NOT dispose of the device with household waste. The Device should be recycled. The Device contains a Lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.

WARRANTY

The device is warranted to be free from defects in material and workmanship for 24 months from its original purchase date under normal use. This warranty extends only to the original retail purchaser with original proof of purchase and only when purchased from an authorized retailer or reseller. This warranty only extends if the device is used in conjunction with authorized Device accessories. If the Device should become defective within the warranty period, contact the retailer that you purchased the device. This warranty does not cover Devices damaged by:

- Accident, misuse, abuse, or alteration.
- Use with unauthorized accessories or use other than as instructed.
- Connecting to incorrect current and voltage.

myBlend shall not be responsible for any incidental, special or consequential damages resulting from the use of this Device. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to two years from date of original purchase. This warranty gives you specific legal rights, and you may also have other rights that vary from country to country. Some countries do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you.

This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose are hereby disclaimed.

If the Device is defective within the warranty period please contact the Customer Service Team. To register your Device and be kept up to date with new and exciting developments from the Customer Service team, please take a moment to sign up on our website www.my-blend.com. To complete your registration online you will need to enter the Device batch number.

SAFETY

SAFETY STANDARDS AND REGULATIONS

The Device has been tested to or complies with the following standards:

- IEC 60601-1:2005+A1:2012+A2:2020. Medical electrical equipment. General requirements for basic safety and essential performance.
- IEC 60601-1-2:2014 Medical electrical equipment. General requirements for basic safety and essential performance. Collateral standard. Electromagnetic compatibility. Requirements and tests.
- IEC 60601-1-6:2010/AMD2:2020 Medical Electrical Equipment – Part 1-6: General Requirements for Basic Safety and Essential Performance – Collateral Standard: Usability.
- IEC/EN 60601-1-11:2015 Medical electrical equipment -- Part 1-11: General requirements for basic safety and essential performance -- Collateral standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment.
- IEC 60601-2-57:2023 Medical electrical equipment. Particular requirements for the basic safety and essential performance of non-laser light source equipment intended for therapeutic, diagnostic, monitoring and cosmetic/aesthetic use.
- IEC 62471:2006 Photobiological Safety of Lamps and Lamp Systems.
- IEC 62366-1:2015+AMD1:2020 Application of usability engineering to medical devices
- ISO 10993-1:2018 Biological evaluation of medical devices - Part 1: Evaluation and testing within a risk management process
- ISO 10993-5:2009 Part 5: Tests for in vitro cytotoxicity.
- ISO 10993-10:2021 Part 10: Tests for skin sensitization.
- ISO 10993-23:2021 Part 23: Biological evaluation of medical devices - Tests for irritation
- ISO 15223-1:2021 Medical devices - Symbols to be used with information to be supplied by the manufacturer - Part 1: General requirements
- ISO 14971:2019 Medical devices - Applications of risk management to medical devices
- IEC/EC 62133-2:2017. Secondary cells and batteries containing alkaline or other non-acid electrolytes - Safety requirements for portable sealed secondary cells, and for batteries made from them, for use in portable applications - Part 2: Lithium systems
- Restriction of Hazardous Substances (RoHS) Directive (2011/65/EU)

- Registration, Evaluation, Authorisation and Restriction of chemicals (REACH) directive (EC 1907/2006)

FAQ

For more information, visit the FAQs section on our website.

SYMBOLS

This symbol marks the placement of the Serial Number

This symbol marks the placement of the LOT Number

This symbol marks the placement of the products manufacturing date

This symbol warns to keep the device dry

This symbol marks that the packaging can be recycled

This symbol means the upper and lower operating temperatures that the device can be safely used at.

This symbol means the upper and lower operating humidity that the device can be safely used at.

This symbol means that the product should be used with caution

This symbol means that anyone who would like to use the product should first read the instructions for use

This symbol marks that the product should not be thrown away with household waste

This symbol marks the manufacturer of the product

The Triman symbol is used to specify which components of the product can be recycled

This symbol signifies that the product adheres to the Restriction of Hazardous Substances directive.

The CE mark indicates that a product complies with European Union (EU) safety, health, and environmental protection requirements

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- IEC 60601-1:2005+A1:2012+A2:2020. Medical electrical equipment. General requirements for basic safety and essential performance.
- IEC 60601-1-2:2014 Medical electrical equipment. General requirements for basic safety and essential performance. Collateral standard. Electromagnetic compatibility. Requirements and tests.
- IEC 60601-1-6:2010/AMD2:2020 Medical Electrical Equipment – Part 1-6: General Requirements for Basic Safety and Essential Performance – Collateral Standard: Usability.
- IEC/EN 60601-1-11:2015 Medical electrical equipment -- Part 1-11: General requirements for basic safety and essential performance -- Collateral standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment.
- IEC60601-2-57:2023 Medical electrical equipment. Particular requirements for the basic

safety and essential performance of non-laser light source equipment intended for therapeutic, diagnostic, monitoring and cosmetic/aesthetic use.

- IEC 62471:2006 Photobiological Safety of Lamps and Lamp Systems.
- IEC 62366-1:2015+AMD1:2020 Application of usability engineering to medical devices
- ISO 10993-1:2018 Biological evaluation of medical devices - Part 1: Evaluation and testing within a risk management process
- ISO 10993-5:2009 Part 5: Tests for in vitro cytotoxicity.
- ISO 10993-10:2021 Part 10: Tests for skin sensitization.
- ISO 10993-23:2021 Part 23: Biological evaluation of medical devices - Tests for irritation
- ISO 15223-1:2021 Medical devices - Symbols to be used with information to be supplied by the manufacturer - Part 1: General requirements
- ISO 14971:2019 Medical devices - Applications of risk management to medical devices
- IEC/EC 62133-2:2017. Secondary cells and batteries containing alkaline or other non-acid electrolytes - Safety requirements for portable sealed secondary cells, and for batteries made from them, for use in portable applications - Part 2: Lithium systems
- Restriction of Hazardous Substances (RoHS) Directive (2011/65/EU)
- Registration, Evaluation, Authorisation and Restriction of chemicals (REACH) directive (EC 1907/2006)

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- IEC 60601-1:2005+A1:2012+A2:2020. Medical electrical equipment. General requirements for basic safety and essential performance.
- IEC 60601-1-2:2014 Medical electrical equipment. General requirements for basic safety and essential performance. Collateral standard. Electromagnetic compatibility. Requirements and tests.
- IEC 60601-1-6:2010/AMD2:2020 Medical Electrical Equipment – Part 1-6: General Requirements for Basic Safety and Essential Performance – Collateral Standard: Usability.
- IEC/EN 60601-1-11:2015 Medical electrical equipment -- Part 1-11: General requirements for basic safety and essential performance -- Collateral standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment.
- IEC 60601-2-57:2023 Medical electrical equipment. Particular requirements for the basic safety and essential performance of non-laser light source equipment intended for therapeutic, diagnostic, monitoring and cosmetic/aesthetic use.
- IEC 62471:2006 Photobiological Safety of Lamps and Lamp Systems.
- IEC 62366-1:2015+AMD1:2020 Application of usability engineering to medical devices
- ISO 10993-1:2018 Biological evaluation of medical devices - Part 1: Evaluation and testing within a risk management process
- ISO 10993-5:2009 Part 5: Tests for in vitro cytotoxicity.
- ISO 10993-10:2021 Part 10: Tests for skin sensitization.
- ISO 10993-23:2021 Part 23: Biological evaluation of medical devices - Tests for irritation
- ISO 15223-1:2021 Medical devices - Symbols to be used with information to be supplied by the manufacturer - Part 1: General requirements
- ISO 14971:2019 Medical devices - Applications of risk management to medical devices
- IEC/EC 62133-2:2017. Secondary cells and batteries containing alkaline or other non-acid

electrolytes - Safety requirements for portable sealed secondary cells, and for batteries made from them, for use in portable applications - Part 2: Lithium systems

- Restriction of Hazardous Substances (RoHS) Directive (2011/65/EU)
- Registration, Evaluation, Authorisation and Restriction of chemicals (REACH) directive (EC 1907/2006)

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Light Tree Ventures Europe B.V.